

Guideline to the PETLINK™ Proposal

26-Mar-2013

Revision J1

PETLINK™ - A proposed
digital interconnect standard
for data acquisition in nuclear
medicine.

A special thanks to the many kind people who have contributed to this document over an estimated 15 year period – 1998 to 2013. These include: L. Byars, J. Reed, P. Luk, E. Breeding, J. Bao, J. Everman, A. Moor, M. Casey, R. Nutt, & B. Bendriem. WFJ

William F. Jones
Senior Key Expert
Siemens Medical Solutions, USA
Molecular Imaging
MED MI ST DPT
810 Innovation Drive
Knoxville TN 37932
voice: 865 218 2254
fax: 865 218 3010
email: william.f.jones@siemens.com
web: www.medical.siemens.com

Table of Contents:

1. Hardware
2. Stream Protocol
3. 64-Bit Packet Format
4. 32-Bit Packet Format
5. 64-bit Tag Packet Format with 56-bit Payload
6. TAXI Adaptations
7. G-Link Adaptations
8. Hardware Extensions for Flow Control

Revision History:

Rev:	Date; Brief Summary of Changes;
G9	30-Jan-2006;
H1	23-May-2007; TAG 3 Formats Expanded; TAG 4 Acquisition Flag; G-LINK History;
H2	1-Nov-2008; Sec. 2 IDLE; Sec 4.2 30-bit Bin-Address Field; TAG 1 Dead-Time; Tag 2 Horizontal Bed Moving Bit; TAG 3 Exp. Format #2;
H3	29-Nov-2010; Acquisition Flag Packets: Time Sync. & Redundant Insertion;
H4	23-Feb-2012; Add "CONFIDENTIAL" to Each Page;
J1	2-Jan-2013; Remove "CONFIDENTIAL"; Sec. 3.1 Filler Packet, Add 2 TOF bits; Sections 3 & 5 Tag_56PL;

Guideline to the PETLINK™ Proposal

26-Mar-2013 Rev J1 Page 1

This document is an evolving guideline for a proposed digital interconnect standard for data acquisition in nuclear medicine - PETLINK™. This proposal is intended primarily for PET and the combination of PET and SPECT applications.

1. Hardware

Current design efforts focus primarily on Fibre Channel (FC) 10-bit transmission characters using FPGA-resident SER/DES (e.g. Xilinx Rocket I/O) and fiber-optic transceivers (e.g. Avago AFBR-57R5AEZ.) [See Fibre Channel web site: <http://www.fibrechannel.org/>]

2. Stream Protocol

Unlike more complex FC protocols, the PETLINK FC stream protocol will emphasize simplicity. No FC frame delimiters will be used. The primary FC ordered set will be the IDLE primitive. IDLE will be used to fill the stream during periods of non-peak throughput. For both 32-bit and 64-bit packet transmission, contiguous transmission characters forming 32-bit and 64-bit data packets will be inserted into this IDLE stream. To help ensure receiver synchronization, the PETLINK transmitter should ensure that IDLE transmission characters always make up some fraction (however small) of the stream - e.g. no fewer than 1 IDLE ordered set for every 128 set-of-4 transmission characters for data even during periods of full-rate data transmission.

3. 64-bit Packet Format

3.1 Following is a 64-bit packet format for detector pair and tag packets:

BIT#	First 32-bit word	BIT#	Second 32-bit word
0-7	AX0-7 (TW 0-7)	0-7	BX0-7 (TW 16-23)
8-15	AY0-7 (TW 8-15)	8-15	BY0-7 (TW 24-31)
16-18	XE0-2 (TW 32-34)	16-18	XE3-5 (TW 44-46)
19-21	AE0-2 (TW 35-37)	19-21	BE0-2 (TW 47-49)
22-24	AI0-2 (TW 38-40)	22-24	BI0-2 (TW 50-52)
25-27	TF0-2 (TW 41-43)	25-27	TF3-5 (TW 53-55)
28	TF6	28	TF7
29	Reserved (Typ. Set to 0.)	29	Reserved (Typ. Set to 0.)
30	Tag_64	30	Prompt / Tag_56PL
31	PS0 = 0	31	PS1 = 1

Where:

AX, BX	Transaxial Head Detector Index (A Head is a portion of a full ring of detectors)
AY, BY	Axial Head Detector Index
XE	Transaxial Encoding (Typically adjacent Heads are not tested for coincidence.)
AE, BE	Energy Window
AI, BI	Depth of Interaction
TF	Time of Flight (2's comp.; 6-bit, TF0-5 or 8-bit, TF0-7; Positive towards "A".)
TW	Tag Word Field: Bits TW 0 – 31 when Tag_64 = 1 & Tag_56PL = 0; Bits TW 0 – 55 when Tag_64 = 1 & Tag_56PL = 1;
Tag_64	Indicates non-event (Tag) 64-bit packet when set to 1; Event 64-bit packet when set to 0. [See also Prompt / Tag_56PL.]
Prompt/Tag_56PL	When Tag_64 = 0, Delayed event if set to 0; Prompt event if set to 1. When Tag_64 = 1, 32-bit Tag packet Pay Load (TW 0-31) if set to 0; 56-bit Tag packet Pay Load (TW 0-55) if set to 1.
PS	Packet Sync; PS0 & PS1.

Note that the two Packet Sync bits are for synchronization with the 32-bit word pair. A receiver of this stream should continuously synchronize on PS0 = 0 in the first 32-bit word. Note also the definition of a 64-bit "Filler" Packet – i.e. a packet to be ignored with "all bits set, save PS0" or "FFFFFFFF7FFFFFFFFF".

A Comment About the 64-bit Filler Packet:

To amplify, this Filler Packet has been defined specifically to aid in the clean termination of packet (data) acquisitions. Typically when a data acquisition is terminated due to elapsed time, such Filler Packets are employed to ensure that any data buffers which may be only partially filled within DMA hardware are properly flushed. This method ensures that whatever (non-filler) packet data had already been received by the elapsed time mark is actually delivered out via DMA into the PC operating system (e.g. Windows) environment - e.g. into a list-mode file. In a typical usage, the user might encounter a few hundred such Filler Packets at the end of a list-mode file. As mentioned, this Filler Packet is intended to be nothing more than a place holder, contains no useful information and may be readily ignored. As of 23-Aug-2012, the intent is for such Filler Packets to only be used with the 64-bit detector-pair packet format described in this Section 3.1.

3.2 Following is a 64-bit packet format for absolute bin address and tag packets:

BIT#	First 32-bit word	BIT#	Second 32-bit word
0-15	BA0-15 (TW 0-15)	0-15	BA16-31 (TW 16-31)
16-19	BA32-35 (TW 32-35)	16-19	BA36-39 (TW 44-47)
20-27	SF0-7 (TW 36-43)	20-27	SF8-15 (TW 48-55)
28-29	Reserved (Typ. Set to 0.)	28-29	Reserved (Typ. Set to 0.)
30	Tag_64	30	Prompt / Tag_56PL
31	PS0 = 0	31	PS1 = 1

Where:

BA	Bin Address Field, Bits 0-39. Projection Space Index.
SF	Scale Factor, Bits 0-15. For, Normalization, Dead-time, etc.
TW	Tag Word Field: Bits TW 0 – 31 when Tag_64 = 1 & Tag_56PL = 0; Bits TW 0 – 55 when Tag_64 = 1 & Tag_56PL = 1;
Tag_64	Indicates non-event (Tag) 64-bit packet when set to 1; event 64-bit packet when set to 0. [See also Prompt / Tag_56PL.]
Prompt / Tag_56PL	When Tag_64 = 0, Delayed event if set to 0; Prompt event if set to 1. When Tag_64 = 1, 32-bit Tag packet Pay Load (TW 0-31) if set to 0; 56-bit Tag packet Pay Load (TW 0-55) if set to 1.
PS	Packet Sync; PS0 & PS1.

Note that the two Packet Sync bits are for synchronization with the 32-bit word pair. A receiver of this stream should continuously synchronize on PS0 = 0 in the first 32-bit word.

4. 32-Bit Packet Format

4.1 Overview

MSB	TXXX XXXX XXXX XXXX XXXX XXXX XXXX XXXX	LSB	T = Tag Bit (1-Tag/0-Event) X = Not Restricted
	0XXX XXXX XXXX XXXX XXXX XXXX XXXX XXXX		Event Packet (Not a Tag Packet)
	10XX XXXX XXXX XXXX XXXX XXXX XXXX XXXX		TAG 1: Time Marker/ Dead-Time Tracking
	110X XXXX XXXX XXXX XXXX XXXX XXXX XXXX		TAG 2: Gantry Motions & Positions
	1110 XXXX XXXX XXXX XXXX XXXX XXXX XXXX		TAG 3: Patient Monitoring: (Gating/Physiological /Head Tracking)
	1111 XXXX XXXX XXXX XXXX XXXX XXXX XXXX		TAG 4: Control / Acquisition Parameters

4.2 32-Bit Packet Format Details

Event Packet: (Most Significant Nibble = 0uuu)

PET Bin-Address Event Packet

0PBB BBBB BBBB BBBB BBBB BBBB BBBB BBBB
| 29 | 0 |

P: Prompt (1 - Prompt; 0 - Delay)
B: Bin Address: 0-29.

Note that older PET systems (< 1998?) designated bit 29 (i.e. B29) as a Transmission/Window (W) bit.

SPECT Event Packet (Obsolete)

0uuu uuH XXXX XXXh SYYY YYYY EEEE EEEE
| 6 | 0 | | 6 | 0 | | 7 | 0 |

u: Undefined (Default: Set to zero.)
H: Head Index (A=0;B=1)
X: Transaxial Xtal Index: 0-6
h: Hit bit (always = 1)

S: Shape bit (Decay: Fast =1 / Slow = 0)
Y: Axial Xtal Index: 0-6
E: Energy Index: 0-7

Guideline to the PETLINK™ Proposal

26-Mar-2013 Rev J1 Page 4

Position Profile Data (Maintenance)

0uuu	uHHH	BBBB	BBBh	XXXX	XXXX	YYYY	YYYY
210	6	0	7	0	7	0	

u: Undefined (Default: Set to zero.)

H: Head number: 0-2

B: Block position: 0-6

h: Hit bit (always = 1)

X: X position data: 0-7

Y: Y position data: 0-7

Energy Discrimination Data (Maintenance)

0uuu	uHHH	BBBB	BBBh	FFFF	FFFF	SSSS	SSSS
210	6	0	7	0	7	0	

u: Undefined (Default: Set to zero.)

H: Head number: 0-2

B: Block position: 0-6

h: Hit bit (always = 1)

F: First energy ADC sample: 0-7

S: Second energy ADC sample: 0-7

Crystal Energy Data (Maintenance)

0uuu	uHHH	BBBB	BBBh	EEEE	EEEE	YYYY	XXXX
210	6	0	7	0	3210	3210	

u: Undefined (Default: Set to zero.)

H: Head number: 0-2

B: Block position: 0-6

h: Hit bit (always = 1)

E: Second ADC sample: 0-7

Y: Y position in block: 0-3

X: X position in block: 0-3

Crystal Time Data (Maintenance)

0uuu	uHHH	BBBB	BBBh	TTTT	TTTT	YYYY	XXXX
210	6	0	6	0	3210	3210	

u: Undefined (Default: Set to zero.)

H: Head number: 0-2

B: Block position: 0-6

h: Hit bit (always = 1)

T: TDC value: 0-6

Y: Y position in block: 0-3

X: X position in block: 0-3

4.2.1 TAG 1: Elapsed Time Marker/Dead-Time Tracking (Most Significant Nibble = 10uu)

Elapsed Time Marker

```

100M MMMM MMMM MMMM MMMM MMMM MMMM MMMM
  | 28                                         0 |

```

M: Elapsed milliseconds: 0-28

Dead-Time Tracking (Legacy Block Singles)

```

101B BBBB BBBB BSSS SSSS SSSS SSSS SSSS
  | 9          0 | 18                      0 |

```

B: Block Number : 0-9

S: Singles per second (per 0.25 sec. May-02; per 0.125 sec. Aug-07): 0-18

Note: Originally this 19-bit S field was defined to represent singles per second. As of May 2002, a change was implemented by Siemens/CTI to down shift by 2 bits the singles/sec value prior to loading into this field – making the S value represent Singles/quarter-second. To get a value in Singles/sec, simply up shift the S value by 2 bits – i.e. multiply by 4. This down shift was done in order to extend the top-end range supported by this 19 bit field. In other words, later versions of the blocks began to effectively detect more than 500k Singles/sec. This 2-bit-down shift allowed support for singles/sec rates for each block as high as 2.1 M/sec and with a negligible loss in precision. As of August 2007 and along similar lines, another change was implemented by Siemens for some later version systems to instead down shift by 3 bits. So for those later version systems to get a Singles/sec value from the S field, multiply S by 8.

Guideline to the PETLINK™ Proposal

26-Mar-2013 Rev J1 Page 6

Expanded Format: Dead-Time Tracking

```
101T TTDD DDDD DDDD DDDD DDDD DDDD DDDD
  2 10|25                                0|
```

D: Data Field: 0-25

T: Type field: 0-2

T field value (0 – 7; one of eight) indicates type of D field content.

Guidelines: T = 0 indicates Block Singles content.

T = 1 – 5 reserved.

T = 6 indicates Lost (Prompt or Delayed) Event Count Packets inserted by the “Second Lossy Node” - e.g. the PDR card or the PSB card.

T = 7 indicates Lost (Prompt or Delayed) Event Count Packets inserted by the “First Lossy Node” - e.g. the Coincidence Processor.

For T = 0: Block Singles Type.

Define whole tag packet as:

```
1010 00BB BBBB BSSS SSSS SSSS SSSS SSSS
  2 10|6      0|18                                0|
```

where: S: Singles per second (per 0.25 sec. May-02; per 0.125 sec. Aug-07): 0-18

B: Block Number: 0-6

Note that this Block Number field is reduced to 7 bits compared to the original 10-bit Block Number for the Legacy case above. Since the range of block numbers reported was always small (< 128), this reduction to 7 bits should be backward compatible.

For T = 6 or 7: Lost Event Counter Type.

Define whole tag packet as:

```
101T TTrr rrrr LLLL LLLL LLLL LLLL LLLL
  2 10      |19                                0|
```

where: T: Type Field: 0-2. Values of 6 or 7 only.

r: reserved.

L: Lost Event Packets: 0-19. The L value represents a count of event packets observed to be lost during the processing of 1M (or precisely 1048575) arriving event packets. In other words, 1M packets come in and only 1M – L packets are transmitted. After each arriving 1M packets are processed, the loss count (L) is transmitted in this tag packet format. If L is zero, this indicates there were no lost event packets for the 1M arriving event packets which were just processed.

Example 32-bit Type-6 (T=6) Lost Event Tag packet with L = 0: **B800 0000** (hex)

Example 32-bit Type-7 (T=7) Lost Event Tag packet with L = 0: **BC00 0000** (hex)

With this Expanded Format, an opportunity exists to also incorporate information related to lost event packets and to lay a foundation for future expansion for on-line reporting of dead-time issues. In regard to Type-7 (T = 7) & Type-6 (T = 6) packets, during periods in which data packet rates in the stream exceed the capacity of the fiber-optic modules – i.e. “Lossy Nodes” such as the Coincidence Processor (GIM) or the PDR card, packets can be discarded (i.e. lost) due to a lack of sufficient communication or processing bandwidth. To correct for these losses, these specific tag packets report the number of event packets which were discarded. [Note that these packets lost due to reasons of limited bandwidth are not to be confused with those packets discarded as being outside the chosen PET field of view (FOV). Those FOV-discarded packets are not to be included in these specific lost-packet calculations.] These lost-event tag packets can be inserted to allow down stream processing (whether on-line hardware or list-mode-driven batch software) of the data stream to correct for these losses as part of an over-all dead-time correction. Such tag packets may not be lost during high rate periods if priority is given by each Lossy Node to lose no tag packets at the expense of lost event packets – i.e. always discard event packets before discarding tag packets.

Guideline to the PETLINK™ Proposal

26-Mar-2013 Rev J1 Page 7

4.2.2 TAG 2: Gantry Motions & Positions (Most Significant Nibble = 110u)

Detector Rotation Position

To be extracted and applied by the Rebinner in PET/SPECT PET mode.

```
1100 0000 BAFF FFFF FFFF FFFF PPPP PPPP
      | 13          0 | | 7      0 |
```

P: PET Rotational Position: 0-7 [For PET mode need 0 - 255d for 0-360° - wfj]

F: Full Precision Rot. Position: 0-13 [Useful for SPECT]

A: CW Rotation: 0 [1 = clockwise gantry rotation as viewed from bed side]

B: CCW Rotation: 0 [1 = counterclockwise gantry rotation as viewed from bed side]

Note: (A=0 & B=0) indicates non-rotating gantry.

Head I "A" Radial Position

```
1100 0001 uuuu uuuu uuR RRRR RRRR RRRR
      | 12          0 |
```

u: Undefined (Default: Set to zero.)

R: Radial Position: 0-12

Head II "B" Radial Position

```
1100 0010 uuuu uuuu uuR RRRR RRRR RRRR
      | 12          0 |
```

u: Undefined (Default: Set to zero.)

R: Radial Position: 0-12

Vertical Bed Position

```
1100 0011 uuuu uuuu uuV VVVV VVVV VVVV
      | 13          0 |
```

u: Undefined (Default: Set to zero.)

V: Vertical Position: 0-13

Horizontal Bed Position

```
1100 0100 uuM HHHH HHHH HHHH HHHH HHHH
      | 19          0 |
```

u: Undefined (Default: Set to zero.)

H: Horizontal Position: 0-19 [Note: For VG30 on Nov-2010 and later mCT systems, this 20-bit H value is always a 2's complement number. When H is near zero - e.g. H = -1(dec) or FFFF(hex), the bed approaches (or is) fully retracted from the FOV. When H is at its largest negative extent, the bed is fully extended into the FOV. The LSB of H represents 0.001cm of horizontal bed movement. In typical practice, the actual precision reported by the PHS has been limited to 0.05cm – i.e. the minimum increment/decrement change for H is 50(dec).]

M: Moving. Single Bit. Set to 0 if horizontal bed motion is unchanging. Set to 1 if horizontal bed motion is changing. [This bit was requested but implementation may be pending.]

Example 32-bit Horiz. Bed Pos. Tag packet (hex) with H = 0: **C400 0000** (hex)

Gantry Left/Right Position

```
1100 0101 uuuu uuuu uuP PPPP PPPP PPPP
      | 12          0 |
```

u: Undefined (Default: Set to zero.)

P: Gantry Left/Right Position: 0-12 [-2287 to 508 decimal in 0.1 mm steps]
(Positive to the left as viewed from bed side of gantry)

Guideline to the PETLINK™ Proposal

26-Mar-2013 Rev J1 Page 8

Source Axial Position & Rotation (Point Source, Rotating Rod, etc.)

```
1100 0110 AAAA AAAA AAAA RRRR RRRR RRRR
      |11      0| |11      0|
A: Axial Position: 0-11
R: Rotational Position: 0-11
```

HRRT Single Photon Source Position

```
1100 0111 uuuu HHHH AAAA AAAA RRRR RRRR
      3210 |7      0| |7      0|
u: Undefined (Default: Set to zero.)
H: Head Number: 0-3
A: Axial Position: 0-7
R: Rotational Position (within the Head): 0-7
```

4.2.3 TAG 3: Patient Monitoring: (Gating/Physiological/Head Tracking) **(Most Significant Nibble = 1110)**

Basic Gating/Physiological Format

```
1110 0EEE uuuu uuuu uuuu uuuu uuuu
      2| |0
u: Undefined (Default: Set to zero.)
E: Expansion Format Control: 0-2.
```

Expanded Format #0 (Legacy)

```
1110 0000 rrrr rrrr rrrr rrrr GGGG GGGG
      2| |0              |7      0|
r: Reserved. (Default: Set to zero.)
G: Gating: 0-7
E: Expansion Format Control: 0-2. (The 3-bit E field is set to 0 – as shown.)
```

Here is a refinement of the Gating byte usage – i.e. the set-of-8 G bits shown above:

```
CPDD DDDD
|5      0|
```

where:

- D: Physiological Data (e.g. Respiratory Phase) Field: 0-5
D content only has validity if flag bit P=1.
- P: Physiological (e.g. Respiratory) Flag Bit – High True.
- C: Cardiac “R Wave” Flag Bit – High True. This bit is intended to be the primary trigger control for driving the automatic FPGA-driven Gating Buffer switching in designs such as the Smart DRAM. [Note: Smart DRAM was an alternate hardware architecture for PET data acquisition that was considered (~2003) but one that did not go into production.]

Example Format #0 32-bit Cardiac R-wave Tag packet commonly used on HiRez: **E000 0080** (hex)

Expanded Format #1 (Gating)

```

1110 0001 rrrr rrrr GGGG GGGG GGGG GGGG
 2 |  | 0          | 15          0 |

```

r: Reserved. (Default: Set to zero.)

G: Gating: 0-15

E: Expansion Format Control: 0-2. (The 3-bit E field is set to 1 – as shown.)

This is the proposed usage format for Gating Tag Packet usage – i.e. the set-of-16 G bits shown above:

```

CPPP rrDD DDDD DDDD
 2 |  | 0  | 9          0 |

```

where:

D: Physiological Data Field: 0-9

D field content only has validity if P field content is nonzero.

r: Reserved. (Default: Set to zero.)

P: Physiological Type Field: 0-2

P field value (0 – 7; one of eight) indicates type of D field content.

Guidelines: P = 0 indicates D field has no valid content. (Default value for P.)

P = 2 indicates R Wave Triggers from multiple ECG devices identified by D - e.g. MR-PET.

D = 0 indicates "primary" R Wave Trigger - e.g. MR-supplied ECG.

D = 1 indicates "secondary" R Wave Trigger - e.g. MR-supplied Oxygen Probe.

D = 2 indicates "tertiary" R Wave Trigger - e.g. customer-supplied ECG

P = 1 indicates Respiratory Trigger Content in D.

[For mMR: D=0 for Gate On. D=1 for Gate Off.]

P = 5 indicates Respiratory Phase/Waveform Content in D.

P = 6 indicates Cardiac Phase/Waveform Content in D.

C: Cardiac "R Wave" Flag Bit – High True. This bit is intended to be the primary trigger control for driving the automatic FPGA-driven Gating Buffer switching in designs such as the Smart DRAM. The intent is for the C bit to only be set = 1 when both the P field = 0 and the D field = 0. Also, the intent is for this C bit to only be set = 1 for single ECG environments. For multiple ECG environments - e.g. MR-PET, see P=2 above.

Example Format #1 32-bit Cardiac R-wave Tag packet commonly used on mCT: **E100 8000** (hex)

Expanded Format #2 (Generic Trigger)

```

1110 0010 rrrr rrrr TTTT TTTT TTTT TTTT
 2|  |0          |15          0|

```

r: Reserved. (Default: Set to zero.)

T: Generic Trigger Field: 0-15

E: Expansion Format Control: 0-2. (The 3-bit E field is set to 2 – as shown.)

The purpose of this proposed Generic Trigger tag-packet format is to support generic requirements for stream marking. These Generic Trigger packets are intended to be used to mark the stream for generic events or with generic data. For example:

1. the patient pushes a button (once or many times) such as to quickly generate an inserted tag packet with the 16-bit T15:0 field = 0;
- and/or
2. the technician pushes a button (once or many times) such as to quickly generate an inserted tag packet with the 16-bit T15:0 field = 1;
- etc.

Note that this Generic Trigger tag packet should not be confused with the TAG 4 Acquisition Flag packet.

Expanded Formats #3 - #6 (Reserved)**Expanded Format #7 (Special Research)**

```

1110 0111 rrrr rrrr SSSS SSSS SSSS SSSS
 2|  |0          |15          0|

```

r: Reserved. (Default: Set to zero.)

S: Special Research Field: 0-15

E: Expansion Format Control: 0-2. (The 3-bit E field is set to 7 – as shown.)

The purpose of this proposed Special Research tag-packet format is to support “special” requirements for research-oriented stream marking. These Special Research packets are intended to be used to mark the stream for special events or with unusual data. For example:

1. patient blood pressure monitoring generates an inserted tag packet with the two-bit S15:14 field set to 1 and the 14-bit S13:0 field set to represent the current BP reading;
- and/or
2. patient temperature monitoring generates an inserted tag packet with the two-bit S15:14 field set to 2 and the fourteen-bit S13:0 field set to represent the current temperature reading.
- etc.

The intent is to allow individual PET researchers the freedom and opportunity to employ this specific-yet-flexible tag-packet format – i.e. a format for marking the event stream with as-yet-unidentified time-critical information. Note that this Special Research tag packet should not be confused with the TAG 4 Acquisition Flag packet.

Motion Tracking (Six Tools with Six Degrees of Motion Tracking - Northern Digital Passive Polaris)

1110 1TTT DDDV VVVV VVVV VVVV VVVV VVVV
 210 210 | 20 0 |

T: Tool_#:0-2 (One of six platforms or tools tracked for position)

D: Degree: 0-2 (Designate one of eight reported position values -

D: 0 1 2 3 4 5 6 7
 V: Q0 Qx Qy Qz Tx Ty Tz Erms)

V: Value: 0-20 (Position values - two's complement - See Polaris Manuals)

4.2.4 TAG 4: Control & Acquisition Parameters (Most Significant Nibble = 1111)

Control tag packets are intended for system component control of devices processing the list mode stream. Currently used by the ECAT Rebinner to control or determine local board reset, mode select, and pass all functions. (See ECAT Rebinner User Guide. Most significant byte = 1111 1100.)

Acquisition parameter tag packets are planned. Acquisition parameter tag packets are expected to be located primarily (but not necessarily) at the start of the list mode stream. Example parameters may include date and time of acquisition start. For absolute bin address streams, the acquisition parameters may include projection space details like span, ring difference, sinogram dimensions, sinogram bin sampling size, etc.

TAG 4: Basic Acquisition Flag Format

1111 1111 CCCC CCCC IIII IIII IIII IIII
 | 31 0 |

- 1111 Bits 31- 28. Most Significant Nibble is Header defined by PETLINK to indicate TAG 4 packet.
- 1111 Bits 27 - 24. Second Most Significant Nibble is the Type field for TAG 4 packet. Set to F(hex) for "Flag".
- C Bits 23 - 16. Check sum. Sum of 3 other bytes in this 32-bit packet. The inserter of the packet is responsible for calculating and including the content of this check sum field. The reader of this packet is responsible for confirming that the check sum value is consistent with the whole 32-bit packet before accepting this packet as valid. Here are four example 32-bit Flag Tag packets in hex with correct check sum:
 TAG 4 (F); Type (F); CS (FF) [FF+00+00=FF]; ID (0000); FFFF0000
 TAG 4 (F); Type (F); CS (00) [FF+00+01=100 or 00]; ID (0001); FF000001
 TAG 4 (F); Type (F); CS (01) [FF+00+02=101 or 01]; ID (0002); FF010002
 TAG 4 (F); Type (F); CS (45) [FF+12+34=145 or 45]; ID (1234); FF451234
- I Bits 15 – 0. Flag Identification Field.

Note that the usual techniques may be employed to contain this 32-bit (TW) tag field within a 64-bit packet. Here are some rough suggestions and comments about the use of the I field:

Incrementing Acquisition Flag Packets:

In this approach the 16-bit I field is simply incremented with each subsequent (non-redundant) insertion of an Acquisition Flag packet into the PETLINK stream. For example, the start of the acquisition for the first bed in a multiple bed position study (not continuous bed motion) would have inserted an Acquisition Flag packet with the I field set to 0. Just prior to ending this first bed position acquisition and prior moving the bed to the next bed position, an Acquisition Flag packet with the I field set to 1 is inserted. For the second bed position, this is marked at the beginning with a Flag insertion which has the I field set to 2. This packet is inserted just as this second bed position becomes stationary. Respectively, a Flag packet is inserted with the I field set to 3 at the end of the second bed position acquisition – just prior to moving the bed along to the next stationary position. And so on. [Now (Dec-2009) we begin to restrict this I-field definition. For incrementing and other general use, only the 15 LSB (I14-I0) in the 16-bit I field are allowed such free use with I15 remaining at zero. See following section on time synchronization.]

Time Synchronization via Acquisition Flag Packets:

Requirements have surfaced for time synchronization between PET and MR, for example. For this, the MSB in this 16-bit I field (I15 bit) is determined to be a “Modality” bit. When I15 is cleared to zero, the tag packet is considered as having originated from the PET system – which supports legacy usage. When the Modality, I15 bit, is set to one, the tag packet is considered to have originated from non-PET signaling. For a synchronizing tag packet resulting from a signal from the MR system, the 16-bit I field is to be 1000 xxxx xxxx xxxx - with “x” meaning that the bit state is unspecified here and remains available for user applications. In other words, we designate the 4 MSB (I15-12) in the I field – when set to 1000 (binary) – as designating a time synchronizing tag packet relating to the timing system within the MR companion to PET. This 4 MSB of the I field can be further specified for other applications as needed.

Redundant Insertion of Acquisition Flag Packets:

Good planning dictates that these Acquisition Flag packets should provide for redundancy. Because the fiber-optic Fibre-Channel communication is known to have bit error rates on the scale of 10^{-12} , there is a small but finite chance that inserted packets may be garbled or lost. For the case of coincidence event packets, losing 1 packet out of trillions is hardly significant. By contrast, losing one unique Flag packet in a list-mode stream can mean a lost (if rare) imaging opportunity.

We plan to deal with this possibility by optionally inserting 1, 2 or 4 (identical, redundant and contiguous) packets – but only when the packets are of a crucial nature such as may be the case for these Flag packets. [Note that this non-zero bit error also is the reason for the check sum field as defined above.] In addition, there should be planning in place such that the application software (and even the FPGA hardware) can support multiple, redundant packets for the Flag packet – both at the time of insertion and at the time of responding to stream content. Insertion of multiple, redundant, identical Flag packets into the stream by FPGA hardware should be automatic - resulting in the specified quantity of contiguous Flag packets inserted with no other packet types intervening. From the beginning, application software (or FPGA hardware) which is written (instantiated) to read and respond to this Flag packet content should respond only to the first arrival of a valid packet – i.e. a packet which shows a valid check sum. If other redundant (and valid) identical packets are received contiguously and just after this first valid packet, these subsequent packets (along with all invalid packets) should be ignored.

In other words, if we find we have problems with corrupted packets, we can choose to have multiple, identical, redundant copies of flag packets (with check sum for validation) which are inserted as a contiguous group to mark each intended point of interest in the packet stream. With this redundancy we gain robustness to identify these points of interest even if a few packets are lost or corrupted.

Further discussion (Dec-2009) on the need for such packet redundancy leads to the following comment. Robustness can also be achieved in some cases without requiring multiple, identical, contiguous insertion of flag packets. As an alternate, the system which originally inserts the packet (e.g. GIM) can instead be required to record (in a remotely accessible register) the millisecond elapsed time tick present when this flag packet was inserted. The acquisition system (e.g. ACS) can then be required to read (presumably as the acquisition ends) and store this time offset value for general reference. The thinking is that if such a critical flag packet is lost (and rarely so, it is expected) from the list-mode stream for some reason, an estimate as to its location within the list-mode file can be inferred from knowing the approximate location to within 1 millisecond. Certainly this alternate approach implies that very few (e.g. as in only one) such critical flag packets are inserted per acquisition.

5. 64-bit Tag Packet Format with 56-bit Payload

The purpose of this section is to describe a newer, 64-bit tag packet format – i.e. one which enables a 56-bit payload field (TW 0 – 55) within 64-bit tag packets. This 56-bit payload field is an expansion on the original 32-bit payload limit (TW 0 – 31) as had been the practice with earlier systems.

5.1 64-Bit Tag Packet Format with 56-bit Payload Field: Overview.

```

PPPP PPPP PPPP PPPP PPPP PPPP PPPP PPPP PPPP PPPP PPPP PPPP PPPP
| 55                                                                 0 |

```

Where:

P Expanded 56-bit Payload Field for Tag Word (TW) Bits: 0 – 55

In practice, this 56-bit field is to be subdivided as follows:

```

TTTT TTTT DDDD DDDD DDDD DDDD DDDD DDDD DDDD DDDD DDDD DDDD DDDD
| 55   48 | | 47                                                                 0 |

```

Where:

T Type: An 8-bit Tag Packet Type Field within Bits 48 – 55.
D Data: A 48-bit Tag Packet Data Field within Bits 0 – 47.

5.2 Block Singles 56-bit Payload Field.

As of 14-Feb-2012, the only 64-bit tag packet which has been defined for using the new 56-bit payload field is the block singles tag packet.

```

TTTT TTTT BBBB BBBB BBBB BBBB SSSS SSSS SSSS SSSS SSSS SSSS SSSS SSSS
| 7    0 | | 15                                                                 0 | | 31                                                                 0 |

```

Where:

T Type: T = 0;
B Block Number: 0-15
S Singles per second: 0-31

In general and especially with recent moves to support continuous bed motion studies, a more wide spread use of the 64-bit bin-address and tag packet format is expected. In such instances, this new-style 56-bit payload tag packet may be supported without compromise. However, note that in those instances in which each of these new-style 56-bit payload tag packets for block singles must still be truncated into the old-style 32-bit payload tag packet, the bit loss that is necessary may not allow all bit-field information to be preserved. Such truncation may be required, for example, during the rebinning step when the 32-bit bin-address event packet is still generated from the 64-bit detector-pair event packet format. In such an instance of forced truncation for the B field above, the current plan is that the 10-bit LSB is simply extracted from the above 16-bit B field and is used to fill the 10-bit B field in the old-style 32-bit block singles tag packet. In addition for the case of the S field above, the appropriate portion of the 32-bit S field above is extracted and is used to fill the 19-bit S field in the old-style 32-bit block singles tag packet. With the latest Aug-2007 convention for a 3-bit down shift in the S field, this means that this specific 19-bit portion (bits 3 – 21) of the 32-bit S field above is extracted and used to fill the 19-bit S field in the old-style 32-bit block singles tag packet.

6. TAXI Adaptations

Not all PET applications require high speed or large event packet bit count - e.g. ECAT EXACT, HR+. Also, rotating PET systems can restrict the event stream bandwidth available when crossing the rotational boundary - e.g. infrared slip ring of the ART, 40 foot coax of the E.CAM pulley assembly as used by PET/SPECT. Through the 1990's the AMD serial communication chip family, TAXIchip, was used at CTI PET Systems, Inc. for event stream communications. These implementations were always with 4 byte (32-bit) event packets which did not support (serial) tag packets. Now for PET/SPECT scanner types, a 5 byte (40 bit) packet is detailed - see 6.3 below.

6.1 EXACT (921) 32-Bit TAXI Event Packet Format

1 st Byte		2 nd Byte		3 rd Byte		4 th Byte	
0-5	MP0-5	0	Prompt	0-4	AD0-4	0-4	BD0-4
		1	Delay				
		2	Multiple				
		3	BR0				
		4	AR0				
		5	BD5	5	AP2	5	BP2
6	BR1	6	AD5	6	AP0	6	BP0
7	AR1	7	Scatter	7	AP1	7	BP1

MP refers to the module pair encoding value.

AD, BD refer to "A" and "B" crystal indexes in the transaxial direction.

AP, BP refer to "A" and "B" crystal-ring ("plane") indexes within the block-pair in the axial direction.

AR, BR refer to "A" and "B" block-ring-to-block-ring indexes in the axial direction.

Note: All successive 4 byte packets are separated by at least one TAXI Command Strobe

6.2 EXACT (922, HR, HR+, ART) 32-Bit TAXI Event Packet Format

1st Byte		2nd Byte		3rd Byte		4th Byte	
0-4	AD0-4	0-4	BD0-4	0	Prompt	0-5	MP0-5
				1-3	AR0-2		
				4-6	BR0-2	6	AD5
5-7	AP0-2	5-7	BP0-2	7	Scatter	7	BD5

MP refers to the module pair encoding value.

AD, BD refer to “A” and “B” crystal indexes in the transaxial direction.

AP, BP refer to “A” and “B” crystal-ring (“plane”) indexes within the block-pair in the axial direction.

AR, BR refer to “A” and “B” block-ring-to-block-ring indexes in the axial direction.

Note: All successive 4 byte packets are separated by at least one TAXI Command Strobe

6.3 PET/SPECT (& LSO-ART) 40-Bit TAXI Event Packet Format

1st Byte		2nd Byte		3rd Byte		4th Byte		5th Byte	
0-7	AX0-7	0-7	AY0-7	0-7	BX0-7	0-7	BY0-7	0-1	AE0-1
	(TW0-7)		(TW8-15)		(TW16-23)		(TW24-31)	2-3	BE0-1
								4	AI0
								5	BI0
								6	Prompt
								7	Tag_40

AX, BX refer to the “A” & “B” crystal index in the transaxial direction.

AY, BY refer to the “A” & “B” crystal index in the axial direction.

AE, BE refer to the “A” & “B” detected photon energy levels.

AI, BI refer to the “A” & “B” “depth of interaction” detector index in the radial direction.

Note: Not all successive 5 byte packets need to be separated by a Command Strobe. Here up to 20 packets of 5 bytes may be adjacent in the TAXI stream - i.e. no more than 20 packets may be transmitted without an inserted TAXI Command Strobe.

6.4 Adapting TAXI to FC – i.e. VME to PC

An investment was made to allow an upgrade for older, VMEbus & TAXIchip compatible tomographs. The goal was to enable the use of the newer PC-based FC systems for limited data acquisition when servicing the older tomograph. This 64-bit detector-pair packet format represents one derived to support data incoming from these older TAXI formats into the newer FC-compatible format. [Two examples of custom PC hardware designed to make this TAXI-to-FC adaptation work are the TFA card and the PDT card - that is, only one such card need be added to the PETLINK & FC-enabled PC, not both.]:

BIT#	First 32-bit word	BIT#	Second 32-bit word
0-7	AD0-7 (TW 0-7)	0-7	BD0-7 (TW 16-23)
8-10	AP0-2 (TW 8-10)	8-10	BP0-2 (TW 24-26)
11-15	AR0-4 (TW 11-15)	11-15	BR0-4 (TW 27-31)
16-18	MP0-2	16-18	MP3-5
19-21	AE0-2	19-21	BE0-2
22-24	AI0-2	22-24	BI0-2
25-29	Reserved (Typ. Set to 0.)	25-29	Reserved (Typ. Set to 0.)
30	Tag_64	30	Prompt
31	PS0 = 0	31	PS1 = 1

Where:

AD, BD	Transaxial Detector Indexes. [Only bits 0-5 are used within each 8-bit field.]
AP, BP	Axial Intra-Block-Pair Crystal Ring Indexes.
AR, BR	Axial Block-Ring Indexes. [Only bits 0-2 are used within each 5-bit field.]
MP	Module Pair Number.
AE, BE	Energy Window
AI, BI	Depth of Interaction
TW	Tag Word Field: Bits TW 0 – 31
Tag_64	Indicates non-event (Tag) 64-bit packet when = 1; event packet when = 0.
Prompt	Prompt event if = 1. Delayed event if = 0.
PS	Packet Sync; PS0 & PS1.

Note that the two Packet Sync bits are for synchronization with the 32-bit word pair. A receiver of this stream should continuously synchronize on PS0 = 0 in the first 32-bit word.

7. G-Link Adaptations

Discussions on the G-Link chip sets from Hewlett-Packard led to the following preliminary formats.

7.1 PET/SPECT 32-Bit G-Link Event Packet Format for Bin Address & Tag (Never Implemented.)

1 st 16 Bit Word	2 nd 16 Bit Word
0-15 BA0-15 (TW0-15)	0-12 BA16-28 (TW16-28)
	13 W (or BA29) (TW29)
	14 Prompt (TW30)
	15 Tag_32 (TW31)

Note: Not all successive double-16 bit word packets need to be separated by Fill Frame. Here up to 50 packets of double 16-bit words may be adjacent in the G-Link stream - i.e. no more than 50 packets may be transmitted without an inserted Fill Frame.

7.2 Rotating Tomograph 48-bit G-Link Event Packet Format (Never Implemented?)

1 st 16 Bit Word	2 nd 16 Bit Word	3 rd 16 Bit Word
0-7 AX0-7 (TW0-7)	0-7 BX0-7 (TW16-23)	0-2 AE0-2
8-15 AY0-7 (TW8-15)	8-15 BY0-7 (TW24-31)	3-5 BE0-2
		6-8 AI0-2
		9-11 BI0-2
		12-13 Reserved
		14 Prompt
		15 Tag_48

Note: Not all successive triple-16 bit word packets need to be separated by Fill Frame. Here up to 33 packets of triple 16-bit words may be adjacent in the G-Link stream - i.e. no more than 33 packets may be transmitted without an inserted Fill Frame.

7.3 PET/SPECT 64-bit G-Link Event & Tag Packet Formats (Never Reached Product Stage)

Detector Pair & Tag

1 st 16 Bit Word	2 nd 16 Bit Word	3 rd 16 Bit Word	4 th 16 Bit Word
0-7 AX0-7 (TW0-7)	0-7 BX0-7 (TW 16-23)	0-2 XE0-2	0-2 XE3-5
8-15 AY0-7 (TW8-15)	8-15 BY0-7 (TW 24-31)	3-5 AE0-2	3-5 BE0-2
		6-8 AI0-2	6-8 BI0-2
		9-14 TF0-5	9-13 Reserved
			14 Prompt
		15 Reserved	15 Tag_64

Bin Address & Tag (Never Implemented)

1 st 16 Bit Word	2 nd 16 Bit Word	3 rd 16 Bit Word	4 th 16 Bit Word
0-15 BA0-15 (TW0-15)	0-15 BA16-31 (TW16-31)	0-5 BA32-37	0-2 AE0-2
			3-5 BE0-2
		6-15 Reserved	6-11 TF0-5
			12-13 Reserved
			14 Prompt
			15 Tag_64

Note: Not all successive quad-16 bit word packets need to be separated by Fill Frame. Here up to 25 packets of quad 16-bit words may be adjacent in the G-Link stream - i.e. no more than 25 packets may be transmitted without an inserted Fill Frame.

7.4 “P39” and “Phoenix” 48-bit G-Link Detector-Pair & Tag Packet Format

1 st 16 Bit Word			2 nd 16 Bit Word			3 rd 16 Bit Word		
0-7	AD0-7	(TW0-7)0-7	BD0-7	(TW16-23)		0-5	XE0-5	
8-10	AP0-2	(TW8-10)	8-10	BP0-2	(TW24-26)	6	AE0 (Scatter)	
11-13	AR0-2	(TW11-13)	11-13	BR0-2	(TW27-29)	7	0	
14-15	0	(TW14-15)	14-15	0	(TW30-31)	8-12	TOF0-4	
						13	Reserved	
						14	Prompt	
						15	Tag_48	

Note: Not all successive triple-16 bit word packets need to be separated by Fill Frame. Here up to 33 packets of triple 16-bit words may be adjacent in the G-Link stream - i.e. no more than 33 packets may be transmitted without an inserted Fill Frame.

7.5 Adapting G-Link to FC

The “Phoenix” ECAT product required a G-Link to FC conversion via the PSA card. This 64-bit detector-pair packet format represents one derived to support data incoming from these G-Link-enabled coincidence processors into the newer FC-compatible format. [Two examples of custom PC hardware designed to make this G-Link -to-FC adaptation work are the PSA card and the PDT card - that is, only one such card need be added to the PETLINK & FC-enabled PC, not both.]:

BIT#	First 32-bit word			BIT#	Second 32-bit word		
0-7	AD0-7	(TW 0-7)		0-7	BD0-7	(TW 16-23)	
8-10	AP0-2	(TW 8-10)		8-10	BP0-2	(TW 24-26)	
11-15	AR0-4	(TW 11-15)		11-15	BR0-4	(TW 27-31)	
16-18	XE0-2			16-18	XE3-5		
19	AE0						
20-29	Reserved	(Typ. Set to 0.)		19-29	Reserved	(Typ. Set to 0.)	
30	Tag_64			30	Prompt		
31	PS0 = 0			31	PS1 = 1		

Where:

AD, BD	Transaxial Detector Indexes.
AP, BP	Axial Intra-Block-Pair Crystal Ring Indexes.
AR, BR	Axial Block-Ring Indexes. [Only bits 0-2 are used within each 5-bit field.]
XE	Transaxial Encoding (Module Pair) Number.
AE0	Scatter Bit
TW	Tag Word Field: Bits TW 0 – 31
Tag_64	Indicates non-event (Tag) 64-bit packet when = 1; event packet when = 0.
Prompt	Prompt event if = 1. Delayed event if = 0.
PS	Packet Sync; PS0 & PS1.

Note that the two Packet Sync bits are for synchronization with the 32-bit word pair. A receiver of this stream should continuously synchronize on PS0 = 0 in the first 32-bit word.

8. Hardware Extensions for Flow Control

Basic PETLINK flow control relies on the proposed Serial FPDP Draft 17.1 from Vita:

<http://www.vita.com/projects/v17d1/wgb/c9-Serial-FPDPv07.pdf>

For coincidence event detection packets in PET, non-systematic loss of a small percentage of the total amount of data during high rate acquisition is a concern but usually not considered a major problem. This limited concern is due to the typically large amounts of data involved - many 100s of millions of packets. However, this PETLINK proposal is built on the concept of including non-coincidence-event (tag) packets into the stream. The partial loss of these tag packets during high rate data acquisition is much more critical. Some tag packet losses can be accommodated - e.g. losing a single millisecond tag packet. Other tag losses are not so easily overlooked - e.g. frame offset control tag packets.

(A note about “Replay” and PETLINK. “Replay” typically refers to post-acquisition processing of list-mode data – i.e. for alternate reframing and/or rebinning. While the primary application for PETLINK is for on-line and real-time data acquisition, PETLINK is effective and desirable for many post-acquisition processing tasks (replay). The most common replay task is converting the collected event stream (list mode) data file into projection space via histogramming. In replay, no data loss is acceptable. Where PETLINK is used to transmit a collected list-mode stream to special hardware, a mechanism for flow control of the stream is desirable to avoid data loss.)

Two levels of extension to PETLINK for flow control are proposed. Level 0 is no flow control. Level 1 adds no additional communication hardware and relies on local FIFO loading intelligence. Level 2 adds a separate low bandwidth upstream communication channel (SUSPEND) for node to node control of flow.

8.1 Flow Control Level 0 (No Active Flow Control)

In PETLINK Flow Control Level 0, a downstream PETLINK node has no mechanism to communicate the SUSPEND condition upstream. Also, no attempt is made to favor one type of stream packet for local FIFO loading over another type in the event of a local FIFO overflow.

8.2 Flow Control Level 1 (Selective FIFO Loading – No Active SUSPEND)


In PETLINK Flow Control Level 1, a downstream PETLINK node has knowledge of potential data loss - e.g. local FIFO near full. This Level 1 stream receiving node actively discards received event packets when necessary in order to minimize the likelihood of tag packet loss. Note that Level 1 compliance is accomplished only by modifying the behavior of the local FIFO loading state machine on the PETLINK stream receive node. Also Level 1 flow control is not intended for replay applications. (Note that for on-line real-time applications, tag packets are expected to be limited to low rates of data transmission - i.e. a few k byte/sec.)

8.3 Flow Control Level 2 (SUSPEND Active – No Selective FIFO Loading)

In PETLINK Flow Control Level 2, a physical upstream connection for SUSPEND is provided from the downstream PETLINK node to the upstream PETLINK node. (See section 7.4 for physical implementation of SUSPEND.) A Level 2 stream receiving node has knowledge of potential data loss (local FIFO near full). This downstream Level 2 node actively communicates SUSPEND to the upstream node. No attempt is made by the downstream Level 2 node to favor one type of stream packet for local FIFO loading over another type in the event of a local FIFO overflow. Level 2 is appropriate for replay and most on-line applications.

8.4 Physical Implementation of SUSPEND Communication.

Envision a three node PETLINK configuration which uses SUSPEND for flow control (Level 2).

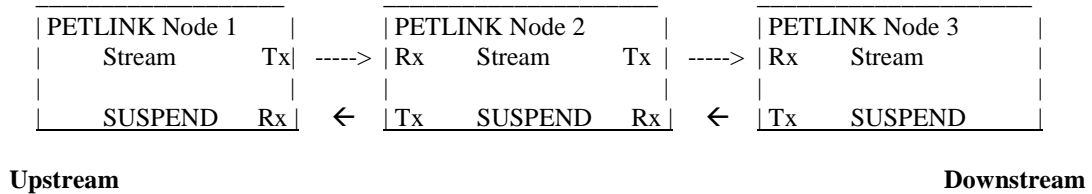

For this proposal, SUSPEND is implemented using the same-type Fibre Channel fiber-optic technology that is used for the downstream communication. The Fibre Channel Ordered Set normally transmitted upstream is the IDLE transmission characters. When a SUSPEND condition is determined by any downstream node, the NOS ordered set (i.e. SUSPEND REQUEST) may be transmitted to all upstream nodes in the chain – i.e. as needed to prevent packet loss.